

Therapy table

Therapy table that is made from 6 mm transparent, impact resistant PETG, with 3 side protectors, and adjustable up to 9 cm in the width.


To be fitted on the armrest upholstery. Adjustable from 28 - 58 mm polstering height with a new designed system.



M1a. Therapy table



M1b. Fitting system

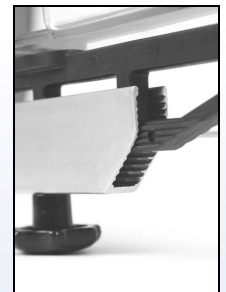

M1c. Detail

Available for 3 widths (internal distance between the armrests):

46-55 cm (seat width of approx. 36-45 cm, depending on wheelchair type)

53-62 cm (seat width of approx. 43-52 cm, depending on wheelchair type)

60-69 cm (seat width of approx. 50-59 cm, depending on wheelchair type)

Therapy table - foldable

The table is produced from impact resistant Polycarbonate and has a frame of chromed steel. The table is fitted with 2 screws under the armrest, and can be adjusted in order to fit all wheelchairs from a seat-width from approx. 35 cm to 50 cm.

When the table is not in use it can be folded along the side of the wheelchair.

The size is 57 x 38 cm (W x D) with a cut-out of 8 cm.



M2a. Therapy table in use



M2b. Therapy table folded