


Standard seat cushions

Standard wheelchair cushion made from high quality cold foam with an especially high load capacity, or with visco-elastic foam. The lower part of the cushion is covered with an anti-slip material. In this way the cushion can be used in any wheelchair without the need to fit with existing velcro strips. The cover is made from Polyester, durable, washable and elastic, and is fitted with a zipper.

The cushion is 5 cm high and all have a seat depth of 43 cm.
Deliverable in 6 widths from 37,5 cm to 50 cm.

Seat cushion with PU foam or high density visco-elastic foam and black polyester or incontinence cover

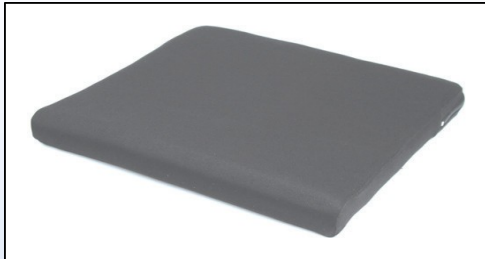

P1. The cushion seen from the front.



P2a. Cushion with incontinence cover



P2b. Anti-slip material

The cushions are available in the respective sizes of both polyurethane foam as high density viscose foam. The standard cushion (5 cm in height), is available with a black polyester or incontinence cover. All models are equipped with zipper and anti-slip bottom.

Note: Available also with individual imprint for an extra charge.

Universal wheelchair cushion (adaptable to the seat depth)



P3a. The cushion with the closed cover



P3b. The cushion with the opened cover and the pre-cutted foam to adapt the seat depth

Universal wheelchair cushion made from high quality foam. The foam comes with precuts so that the cushion is easily adjustable to the depth of the wheelchair seat. The incontinence cover is equipped with velcro stripes for easy closing.

The cushion is available in seat widths from 40 to 51 cm.
The depth is adjustable from 48 cm to 42 cm.

Seat cushion with a double layer (PU foam and visco-elastic foam)



P4a. Incontinence cover



P4b. Anti-slip material



P4c. Viscoelastic foam

Seat cushion made from one layer of normal PU foam and one layer of visco-elastic foam. It comes with an incontinence cover. The cushion is deliverable in 5 different sizes from 40 to 50 cm.



P5. Wheelchair cushion with valve for regulation of air pressure

Mobilex air cushion

The Mobilex air-cushion is equipped with a special valve for air pressure. The user "immerses" in the foam, which conforms to the user and distributes the user's weight over the surface of the cushion. By closing the valve the foam maintains its form for an optimal pressure distribution.

Available in 3 different sizes (W x D x H):

- 40 x 40 x 8 cm
- 45 x 45 x 8 cm
- 50 x 48 x 8 cm



P6. Swivel seat with Ø 45 or 50 cm

Mobilex swivel seat

Can be used in a car or other place where rotability is wanted.

Diameter: 45 or 50 cm.

Max. user weight: 100 kg.



P7. Wedge with incontinence cover

Mobilex wedge with Decutex cover

The wedge cushion is made with 2 layers of foam. The bottom layer is a stable pelletfoam, the upper a viscoelastic foam, which easily adapt.

The Decutex incontinence cover is easy to clean.

Size (W x D x H): 45 X 43 X 10/0 cm



P8. Wedge made of regenerated foam

Mobilex standard wedge

Made of regenerated foam.

Available in 2 sizes (W x D x H):

- 45 X 45 X 8/1 cm
- 45 X 45 X 10/1 cm